

freeway

## Flow Shower Cradle



# Key features

**Head Cushion\***  
Experience additional comfort with the head cushion.

**Harness Points**  
Fixing points enable the attachment of a chest/ sternum harness for even greater security.

**Hand Lever Tilting Mechanism**  
The shower cradle can easily be reclined from 5° to 35°.

**Padded Side Rails\***  
Elevate comfort and support with padded upper and/or lower side rails.

**Height Adjustable Frame**  
The frame can be manually height adjusted to three set heights, in increments of 25 mm.

**Lap Belt\***  
Add an extra layer of security with an optional lap belt.

**Commode Tray And Pan\***  
Easily transform the cradle into a convenient commode.

**Aperture Hole Opening**  
The cradle can have an aperture for use with a commode pan.

**Durable Frame**  
The powder-coated stainless steel frame is strong and stable.

**Anti-Tangle Castors**  
High-quality 125 mm covered braked castors prevent dirt build-up.

**Mesh Fabric**  
The fabric, which allows water to pass through during showering, is removable for easy cleaning and drying.



\*Extra charges apply.

# Weights and dimensions


Cradle size	Small	Medium	Large
SWL	130 kg / (286 lb)	130 kg / (286 lb)	130 kg / (286 lb)
A* – Seat height	592 mm / (23.3 ")	592 mm / (23.3 ")	592 mm / (23.3 ")
B – Cradle end height	176 mm / (6.9 ")	116 mm / (4.6 ")	116 mm / (4.6 ")
C* – Handle height	940 mm / (37 ")	1,050 mm / (41.3 ")	1,050 mm / (41.3 ")
D* – Overall height	1,106 mm / (43.5 ")	1,214 mm / (47.8 ")	1,214 mm / (47.8 ")
E – Base frame length	800 mm / (31.5 ")	800 mm / (31.5 ")	800 mm / (31.5 ")
F – Back length	620 mm / (24.4 ")	750 mm / (29.5 ")	750 mm / (29.5 ")
G – Seat length	300 mm / (11.8 ")	350 mm / (13.8 ")	350 mm / (13.8 ")
H – Leg length	360 mm / (14.2 ")	450 mm / (17.7 ")	450 mm / (17.7 ")
I – End length	150 mm / (5.9 ")	150 mm / (5.9 ")	150 mm / (5.9 ")
J – Internal width	440 mm / (17.3 ")	440 mm / (17.3 ")	520 mm / (20.5 ")
K – Overall width	502 mm / (19.8 ")	502 mm / (19.8 ")	582 mm / (22.9 ")
Upright length	950 mm / (37.5 ")	1,134 mm / (44.6 ")	1,134 mm / (44.6 ")
Reclined length	1,260 mm / (49.6 ")	1,508 mm / (59.4 ")	1,508 mm / (59.4 ")
No. of height settings	3	3	3
Height increments	25 mm / (1 ")	25 mm / (1 ")	25 mm / (1 ")

## Details

Warranty	5 years (excluding fabric)
	2 years (fabric)
Frame material	Powder coated stainless steel
Fabric material	Polyester mesh

## Standards

IEC 62366-1:2015
EN ISO 15223-1:2016
ISO 17966:2016
Fabric material

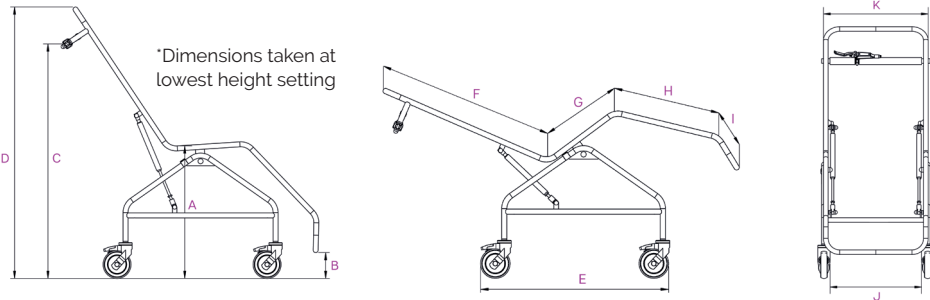

To request a quote or free assessment, please call us on  
**01924 840 100** or visit [www.prismmedical.co.uk/contact-us](http://www.prismmedical.co.uk/contact-us)

**Prism Medical UK Head Office**

Unit 4, Jubilee Business Park,  
Grange Moor, Wakefield  
WF4 4TD

**Prism Medical UK Manufacturing Centre**

Unit 1, Tir Llwyd Industrial Estate,  
Kinmel Bay, Conwy  
LL18 5JZ

Telephone: **01924 840100**

Email: [info@prismmedical.co.uk](mailto:info@prismmedical.co.uk)

Website: [www.prismmedical.co.uk](http://www.prismmedical.co.uk)

Please visit our website to find out more about the  
**Freeway Flow Shower Cradle**

