

Freeway TA180

The Freeway TA180 is a sit-to-stand transfer aid designed for individuals with partial mobility. It enables safe, stable movement over short distances, such as between beds, chairs, and toilets, without the need for a wheelchair.

Key features include padded leg supports, adjustable seat paddles, and locking rear castors, all contributing to user comfort and carer control. The TA180 accommodates a variety of seating configurations and offers interchangeable foot trays to suit individual client needs.

UK Part Number	1301TA-FWTA180 (with standard foot tray) 1301TA107435 (with wide foot tray)
User weight limit (SWL)	185kg / 29st / 407lbs
Total weight (unboxed)	<40kg
Castor size	100mm (rear-braked)
Ground clearance	78mm
Under-bed clearance	128mm
Turning circle	1010mm
Length	937mm
External width	681-1048mm
Internal width	535-888mm
Height	1021mm
Tray width	365mm (standard), 466mm (wide)

Offers optional extra support

an optional support belt clips onto the upper frame for added security during transfer




Promotes user comfort and stability

padded knee support with 50mm cushioning and optional back support belt for added reassurance



Adapts to individual needs

removable foot tray available in standard or wide* options, compatible with various seating setups and supports standing on the floor

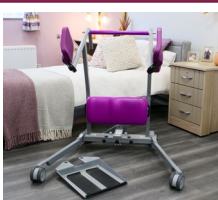

Accommodates chairs with armrests

paddle angle and positioning designed for compatibility with riser-recliner armchairs



Secures the user during transfer

swing-away seat paddles lock into place once upright, reducing fall risk



Functions as a standing frame or rehab tool

can assist with gait training or independent standing when used by trained professionals

Facilitates room-to-room transfers

twin castors (rear braked) allow smooth movement, even over carpeted flooring



* Wide foot tray might not fit standard wheelchairs, 440mm commode chairs, dining chairs, and some riser-recliners.