

# Gentona

## A versatile, variable-height bathing system designed for intensive use in care homes and hospitals.

Designed specifically for high-traffic, wet environments, the Gentona is fully height adjustable and comes with powered transfer seat to deliver exceptional levels of care. With its energy-efficient keyhole shape and a variety of water-saving features, this versatile bathing system is designed to deliver high levels of performance yet it is still extremely economical to run.

Being height adjustable, the Gentona ensures a safe, comfortable working height for care staff, and its easy-to-clean surfaces help to reduce the risk of cross-infection. The powered transfer seat lowers to a convenient height for wheelchair or hoist access and gently moves the bather in and out of the bath.



Lightweight, ergonomic hand operated shower unit

Two position footboard\* customises the bathing experience to the needs of each bather

TMV3 water temperature controls reduces risk of scalding

Lightweight, ergonomic handset with battery back-up system

Bath can be raised or lowered to the optimum working height for carers

Sensory control panel\* Combining hydrotherapy and chromatherapy to enhance bather experience

Folding nursing armrests for enhanced client safety

\*Optional extra


### Options



See pages 6-7.

### Dimensions

CAD drawings available from [www.gainsboroughbaths.com](http://www.gainsboroughbaths.com)

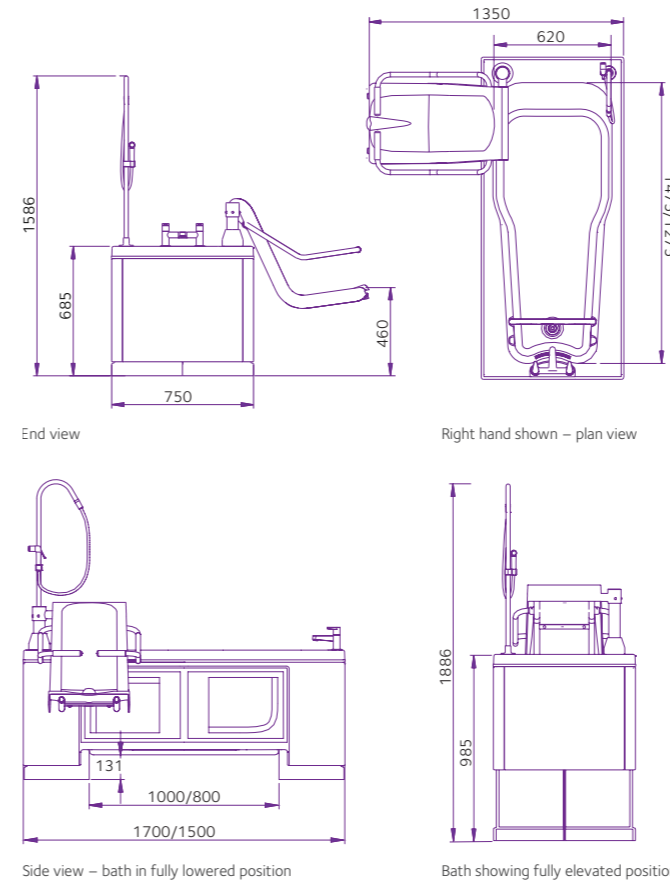

### Key features

- Fully-adjustable, hi-lo bath ensures a safe, comfortable working height for carers
- Powered seat traverse effortlessly transfers the bather in and out of the bath, reducing the need for manual handling
- Innovative keyhole design reduces water usage and helps cut energy costs
- Depth-level indicators allow carers to pre-fill the bath to a required depth, saving water and eliminating unnecessary additional top-ups
- Detachable transfer chair system to assist carers when transporting bathers between rooms (optional)
- The 1700mm full length bath boasts extra-long internal bathing length and wide seating area, ensuring optimum comfort for bathers
- Double skinned bath panels for added durability and strength
- The ultra-compact 1500mm bath is perfect for space restricted environments
- TMV3 controls maintain a safe, comfortable, pre-set water temperature, reducing any risk of scalding
- Stylish, chrome shower system with easy-to-use hand control
- Lightweight, ergonomic handset with battery backup safety system in case of mains power failure
- 160kg (25 stone) safe lifting weight accommodates the majority of bathers
- Left- or right-hand configurations adapt to any bathroom layout



### Care level 6

The Gentona allows full carer assistance with hi-lo facility and powered seat.

